









EXTENDED RANGE



Continuing to expand the Guardian range of single-use products to bring you an exciting array of high quality ancillary items.

Our new products include:

-  **Facemasks** to fulfil all your requirements for surgery and procedures
 -  **Splash Shields** when extra protection is required for your eyes
 -  **Gowns** to keep you dry and comfortable
 -  **Patient Gowns** for modesty and convenience
 -  **Cornea Covers** to prevent corneal drying and abrasion during general anaesthesia
-  *For priority news of the next exciting addition to this range, please register your interest with us for immediate notification: info@bdf.ltd.uk*

In the meantime, please read on ...

Maximum Fluid Protection Anti-fog Facemask



Non-fogging wraparound optical quality eyeshield assures maximum visibility and facial protection from impacting fluids. 2-ply filter, cellulose outer and innerfacing for skin comfort. Comfortable anti-fog foam strip prevents fogging of eyewear. Aluminium wire for secure fit.

(Peach) with Ties & Eye Shield

2-ply filter, cellulose outer and innerfacing, for skin comfort. Comfortable anti-fog foam strip prevents fogging of eyewear. Aluminium wire for secure fit.

(Peach) with Ties

Guardian Code
3505PA

Case Qty
100 (4x25)

Equivalent	
K-Clark	48247; 48248
3M	
J & J	4232
Allegiance	AT74635
NHS	BWM019

Guardian Code
3503PA

Case Qty
300 (6x50)

Equivalent	
K-Clark	48207
3M	
J & J	4234
Allegiance	AT74535
NHS	BWM020

Fluid Protection Anti-fog Facemask

Equivalent	
K-Clark	62114
3M	1818FS
J & J	4232
Allegiance	AT74635; AT74671
NHS	BWM019

Guardian Code
2505BA

Case Qty
100 (4x25)

Non-fogging wraparound optical quality eyeshield assures maximum visibility and facial protection from impacting fluids. An extra protection filter provides excellent resistance to fluid impact. Cellulose outer and innerfacing for skin comfort. Comfortable anti-fog foam strip prevents fogging of eyewear. Aluminium nose wire for secure fit.

(Blue) with Ties & Eye Shield



Guardian Code
2508BA

Case Qty
100 (4x25)

(Blue) with Ear Loops & Eye Shield

Equivalent	
K-Clark	47147
3M	
J & J	4267
Allegiance	AT74631
NHS	

Fluid Protection Anti-fog Facemask



An extra protection filter provides excellent resistance to fluid impact. Cellulose outer and innerfacing for skin comfort. Comfortable anti-fog foam strip prevents fogging of eyewear. Aluminium nose wire for secure fit.

(Blue) with Ties

(Blue) with Ear Loops

Guardian Code
2503BA

Case Qty
300 (6x50)

Equivalent	
K-Clark	62113
3M	1818
J & J	4231; 4234
Allegiance	AT74535; AT74335
NHS	BWM016/020

Guardian Code
2501BA

Case Qty
300 (6x50)

Equivalent	
K-Clark	62115
3M	
J & J	4224; 4229
Allegiance	AT73531
NHS	

All Guardian facemasks shown on this page meet EN 13683:2005 Class II (* except 2510L = Class IIR)



Anti-fog Facemasks

Poly outer and cellulose innerfacing for skin comfort.
Aluminium nose wire for secure fit.

Medical grade tape strip to prevent fogging of eyewear.

Guardian Code

1624G

Case Qty
300 (6x50)

Equivalent	
K-Clark	49215
3M	
J & J	
Allegiance	AT72835
NHS	

Soft foam strip to prevent fogging of eyewear.

(Green) with Ties

Guardian Code

1514G

Case Qty
300 (6x50)

Equivalent	
K-Clark	49214; 49310
3M	
J & J	4231
Allegiance	AT71235
NHS	

EMA film strip to prevent fogging of eyewear

(Green) with Ties

Guardian Code

1504G

Case Qty
300 (6x50)

Equivalent	
K-Clark	49230
3M	
J & J	4238
Allegiance	AT752005
NHS	

Equivalent	
K-Clark	48220
3M	1838
J & J	4233
Allegiance	AT51035
NHS	BWMO26

Guardian Code

1606G

Case Qty
300 (6x50)

Chamber Mask

Duck bill style with large breathing chamber. Poly outer and cellulose innerfacing for skin comfort. A full width nosepiece ensures a proper fit.

(Green) with Ties



Sensitive Skin

Made with no inks or dyes that could potentially irritate the skin. Cellulose outer and innerfacing for skin comfort. The aluminium nose wire gives a secure fit.

(White) with Ties

Equivalent	
K-Clark	49310
3M	
J & J	
Allegiance	AT7305
NHS	

Guardian Code

1502W

Case Qty
300 (6x50)

Equivalent	
K-Clark	47625
3M	1812
J & J	4225
Allegiance	AT75435
NHS	

Guardian Code

2510L

Case Qty
300 (6x50)

Laser Mask * meets EN 13683:2005 Class IIR

With an extra-protection filter rated at 0.1 micron particle size. A cellulose outer and innerfacing for skin comfort. Aluminium nose wire for secure fit.

(Grey) with Ties



All Guardian facemasks shown on this page meet EN 13683:2005 Class II

Equivalent	
K-Clark	48201
3M	
J & J	42280
Allegiance	AT71035
NHS	BWM006

Guardian Code
1503B

Standard Surgical Mask

Cellulose outer and innerfacing for skin comfort

Case Qty
300 (6x50) (Blue) with Ties



Equivalent	
K-Clark	UN48100
3M	1810F
J & J	4230
Allegiance	AT71039
NHS	BWM005

Guardian Code
1403B

Lightweight Surgical Facemask

Our lightest facemask.
Poly outer and cellulose innerfacing for skin comfort..

Case Qty
300 (6x50) (Blue) with Ties



Equivalent	
K-Clark	48390
3M	1818
J & J	
Allegiance	
NHS	BWM016

Guardian Code
1303B Soft Surgical Facemask

Our softest facemask.
Poly outer and bi-component innerfacing for maximum softness.

Case Qty
300 (6x50) (Blue) with Ties



Cone Style Mask

Easy to wear. The non-collapsing, moulded polyethylene breathing chamber allows you to "breathe easy" and gives good protection from air-borne bacteria. The secure closure from the heavy gauge aluminium nose strip prevents leakage.

Guardian Code
1550C

Case Qty
300 (6x50)

Equivalent	
K-Clark	47200
3M	1800+NL
J & J	4235
Allegiance	AT7509
NHS	BWM022

SURGICAL MASKS



All Guardian facemasks shown on this page meet EN 13683:2005 Class II (* except N95 - n/a)

Equivalent	
K-Clark	47080
3M	
J & J	
Allegiance	AT71021

Guardian Code

1001B

Case Qty
300 (6x50)

Standard Procedure Mask

Every-day Procedure Mask Lightweight poly outer and soft, non-irritating tissue innerfacing.

(Blue) with Ear Loops

Equivalent	
K-Clark	62357
3M	
J & J	
Allegiance	AT7511

Guardian Code

1401B

Case Qty
300 (6x50)

Light Procedure Masks

Our lightest mask. Poly outer and cellulose innerfacing for skin comfort.

(Blue) with Ear Loops

Equivalent	
K-Clark	62357
3M	
J & J	
Allegiance	

Guardian Code

1401G

Case Qty
300 (6x50)

(Green) with Ear Loops



Mint Scented Procedure Mask

When you need a "breath of fresh air", this spearmint-scented mask can be a real asset. Poly outer and cellulose innerfacing for skin comfort.

(Mint Green) with Ear Loops

Guardian Code

1401GM

Case Qty
300 (6x50)

Equivalent	
K-Clark	
3M	
J & J	
Allegiance	



Isolation Mask

Lightweight poly outer and soft tissue innerfacing for comfort

(Lemon) with Ear Loops

Guardian Code

1001Y

Case Qty
300 (6x50)

Equivalent	
K-Clark	47117
3M	
J & J	
Allegiance	AT7311-Y

SplashShield

Guardian
Code
52.00.01

Maximum face protection for the wearer

Case Qty
288

- Anti-fog
- Optically clear Impenetrable foam barrier for wide coverage
- A soft strap for comfort and lightweight construction
- Disposable
- Latex free
- Cost efficient
- CE Marked

Guardian
Code
52.00.11

Case Qty
96



The SplashShield comprises an optically clear face screen, with an antifog coating, supported from a foam pad held in place on a users forehead by an elastometric band attached at each end to the optically clear screen adjacent to the foam pad.

The SplashShield is, when flat, about 33cm (13 inches) wide by 19cm (7 ½ inches) deep.

When worn by a user, it covers the face from forehead to chin and, looking from the side, from between the ear and eye on one side to between the ear and eye on the other side.

It has been designed to protect clinicians and other personnel who are at risk from being struck on the face by blood and other human tissue which may potentially be contaminated by untreatable pathogens like HIV, Hepatitis C or prions associated with new variant CJD.

The device has been designed to meet the clarity, splash protection and impact protection requirements of EN 166.

THE FABRIC

Ahlstrom's BVB is a triple-layered fabric, the only one of its kind anywhere, and employs an active fabric technology that features a responsive middle membrane sandwiched between inner and outer fabrics made of continuous fine filaments. The inner layer provides a soft touch to the wearer's skin while the outer layer provides additional repellency and strength. The BVB fabric passes ASTM F1671 and provides the maximum viral protection due to monolithic membrane technology.

THE FABRIC BREATHES

The most significant advantage of this breakthrough fabric is that this monolithic inner membrane automatically responds to the rising body temperature of the wearer, increasing the rate of moisture transfer, helping to keep you comfortable without compromising the viral barrier.

MAJOR FEATURES

- Absolute microbe barrier
- Selective permeability
- High pressure water resistance
- Good tear strength
- High moisture transmission rate
- Water and liquid proof

(The alternative technology offered by other manufacturers uses a micro-porus film, which stretches to create tiny holes that could also allow other molecules to pass through)

BVB provides the best combination of properties for protection and comfort available today.

For detailed information please visit Ahlstrom's web site:

www.viralbarrier.com



THE GOWN


Always proactive in bringing the best technological advances to its Guardian brand, BDF Ltd have taken this remarkable fabric and combined it with our design expertise to create this truly outstanding surgical gown.

Exceeding the stringent demands of the healthcare sector and listening to your concerns, we have enhanced our latest design by welding the sleeve seams rather than sewing - this gives added barrier protection.

Don't compromise - Wearing the new

BREATHABLE, IMPERVIOUS GOWN

from

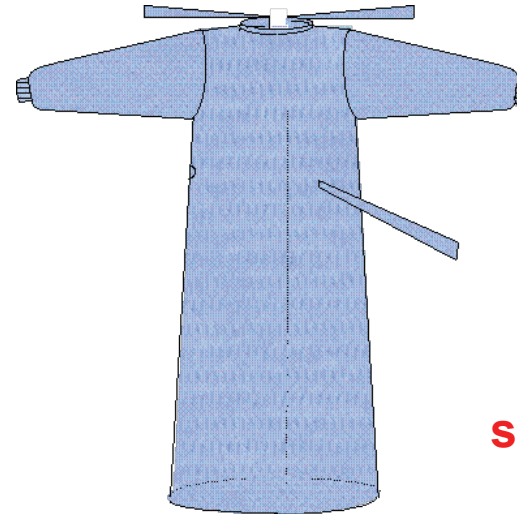

means that you don't have to choose between safety and comfort

Code	Size	Case Qty
11.60.02	Medium	24
11.60.03	Large	24
11.60.04	Extra Large	24

Feels comfy & cool

All-in-one styling, breathable, lightweight, water repellent material

Protects your clothes

Durable plastic film provides zoned reinforcement - just where you need it!

2 sensible sizes

Generously sized, both the Large and Extra Large gowns have a perforated section at the hemline to allow quick and easy adjustment of the length

Easy to dispose of

Made from 100% recyclable material

Guardian Code	Size	Case Qty
51.60.01	Large	80
51.60.02	X-Large	80

STAY DRY!

BE COMFORTABLE

When you wear the

New  *guardian*
Wash Room Gown



Too long?

simply
tear
along
the
perforations



Patient Examination Gown

Ease of use

For both examiner and patient

Comfortable fit

A generously cut, short sleeved, round necked style

Modesty

The discreet, dark colour and conveniently placed neck and waist ties help to create a garment that aids patient composure

The ideal single use patient examination gown

Guardian Code	Size	Case Qty
51.60.03	One-size	100



The Cornea Care eye protection dressing is a clean and effective means of closing the eyes. Comparative mobility of the upper eyelid has been utilised to shut against the lower lid. It can be used on patients undergoing general anaesthesia, long term ventilation etc.

For use by anaesthetist/anaesthesiologist or intensivist to avoid corneal drying and reduce corneal abrasion during general anaesthesia, sedation or similar procedures where corneal drying or injury is a likely concern.

The closed status of the eyes is visible through the clear film, which is made of transparent material, this enables continuous monitoring of the eyes by the user.


Cornea Care aids the formation of a micro-moisture environment underneath the clear film, resulting from the difference between the temperature of the body and the operating theatre. This moisture keeps the eyelid margins moist.

Complete closure of the eyes aids natural lubrication.

Ordinary primitive tapes are porous and exposed to cross contamination from hands, saliva from intubation and anaesthetic drug traces. This is addressed by a no touch technique by providing split backing release paper.

Decrease the loss of eyebrow and eyelash hair upon removal.


Improve the quality of patient care.

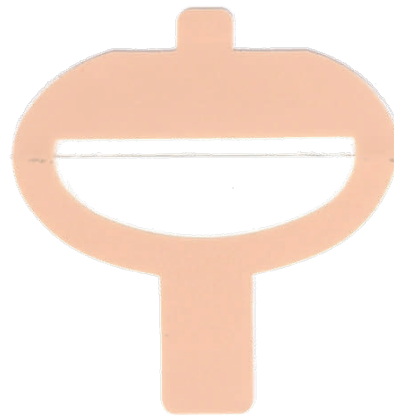


Cornea Care Eye Covers

Guardian
Code
52.00.02
Case Qty
500 prs


CORNEA-CARE instructions for use

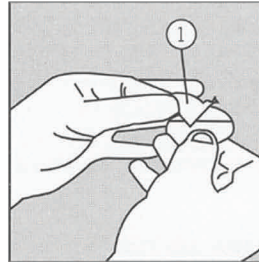


Remove the upper release paper (1) and hold the device in between thumb and index finger.

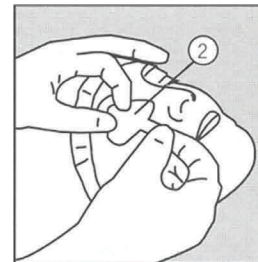

Gently lift eyebrow with an index finger and place the upper adhesive portion (2) over the upper eyelid.

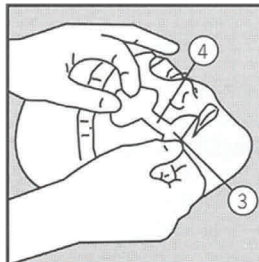

Peel lower release paper (3) and pull tab (4) to apply lower adhesive limb portion over the maxillary area of face.

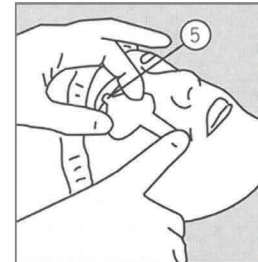

Gently smooth the CORNEA-CARE to obtain eye closure. At the conclusion of the procedure simply pull downward from the upper tab (5).

PRECAUTIONS

Ensure that adhesive surface is **NOT** in contact with cornea. Avoid usage on a patient allergic to tapes/plasters; on steroid therapy or with skin disorders. Additional precautions should be taken when there is a risk of damage to the eyes eg: prone position.



We are convinced that you will be as excited as we are by this innovative product and would be delighted to discuss the many benefits of using this product in greater detail with potential users. Please contact us directly or through your Account Manager, or visit www.bdf.ltd.uk for more information.

cornea-care*

by  AQZ™



what
next ?

WATCH THIS SPACE!

guardian WATCH THIS SPACE!

